

HUNTLEIGH

Baby Dopplex 4000 Series Fetal Monitors User Guide

Refer to instructions for use for detailed set-up, operating, cleaning & maintenance information

For System / Patient Set-up – see reverse of card

Reset Contractions trace baseline



1. Check belt tension
2. Press Contractions Zero key on control panel

Trace annotation

1. To print an annotation message:
 - a. Use Menu key (left of screen) to scroll through message list
 - b. When desired message displayed, press Clinical Event Marker key (on control panel)
2. To print marker without text, press Clinical Event Marker key



Annotation Messages


1. U/S transducer moved
2. Toco transducer moved
3. Mother moved
4. Treatment
5. Vaginal Exam
6. Oxygen
7. Syntocinon
8. Pethidine
9. Epidural
10. Epidural top-up

Clear an alarm

1. Confirm cause & check fetal condition
2. Press soft-key indicated on screen to reset alarm
3. After reset, alarm is still active – if condition persists, alarm will re-trigger after user selected delay
4. To change alarm modes / settings refer to 'Change System Settings'

Fast-feed paper

1. Switch printer off - automatically performs short fast-feed
2. For additional feed with printer off:
Press centre soft-key

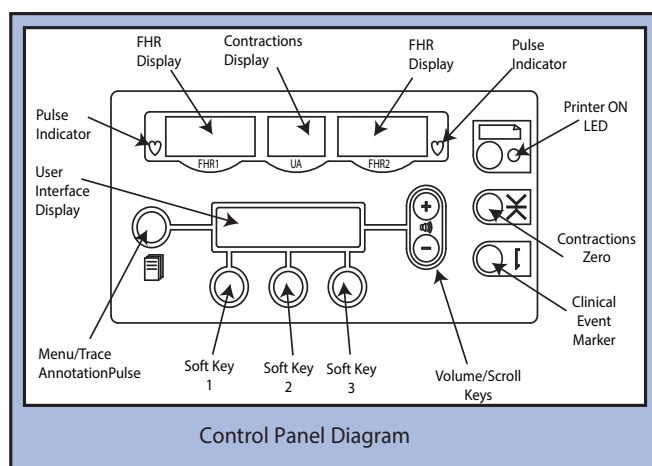

Instrument Set-up

1. Connect power lead (socket on rear panel) & switch unit on (switch above power socket)
2. Check the following: – refer to instructions for use to adjust
 - a. Displayed date & time correct
 - b. Sufficient paper in tray – push tray in firmly on both sides to latch closed
 - c. Chart speed setting correct (shown on display)
 - d. To change system settings see 'Change system settings' below
3. Connect the transducers – note plug/socket colour coding.

Patient / trace Set-up

1. Position transducers & secure with belts
 - a. Ensure sufficient & correct type of gel is used (U/S transducer ONLY).
 - b. Adjust U/S transducer position for clear fetal heart sounds - avoid using umbilical artery signal for best pick-up
 - c. TWINS monitoring – use left & right soft keys under the screen to select FHR1 or FHR2 audio
 - d. Press the 'Contractions Zero' key – sets base line to +20
2. To start or stop printing, press the Printer On/Off key.
The yellow LED lights when the printer is running.



Change system settings

1. Switch printer OFF
2. Press Menu key (to left of screen)
3. Use Volume Up/Down keys to select function
4. Change function settings using soft-keys below screen as indicated
5. Press Menu key to exit
6. Use soft-keys to save changes or quit as indicated

Trouble-shooting

1. Poor FHR trace
 - a. Adjust U/S transducer position for clear heart sounds – avoid using umbilical artery blood flow signal.
 - b. Check sufficient gel present on U/S transducer face
2. Poor contractions trace
 - a. If baseline is off the bottom of the chart, check belt tension & re-zero
 - b. If poor registration, adjust transducer position
3. Refer to instructions for use for further fault-finding assistance.